

Creative

CREATIVE LIFESCIENCES

啟新生技

唯有啟新 · 專業放心

Thermo Scientific
Culture Media in Bioprocess Containers

Safe & simple validation

Thermo
SCIENTIFIC


Where validation meets innovation

Media fills play a critical role in the validation of aseptic manufacturing processes.


To ensure a positive outcome, it is essential to use a sterile growth medium as the placebo in place of your final product.

Stay in total control of this demanding process using compliant, validated Thermo Scientific™ Culture Media.

Supplied, ready-to-use, in a range of convenient, plug-and-play bioprocess containers (BPCs).

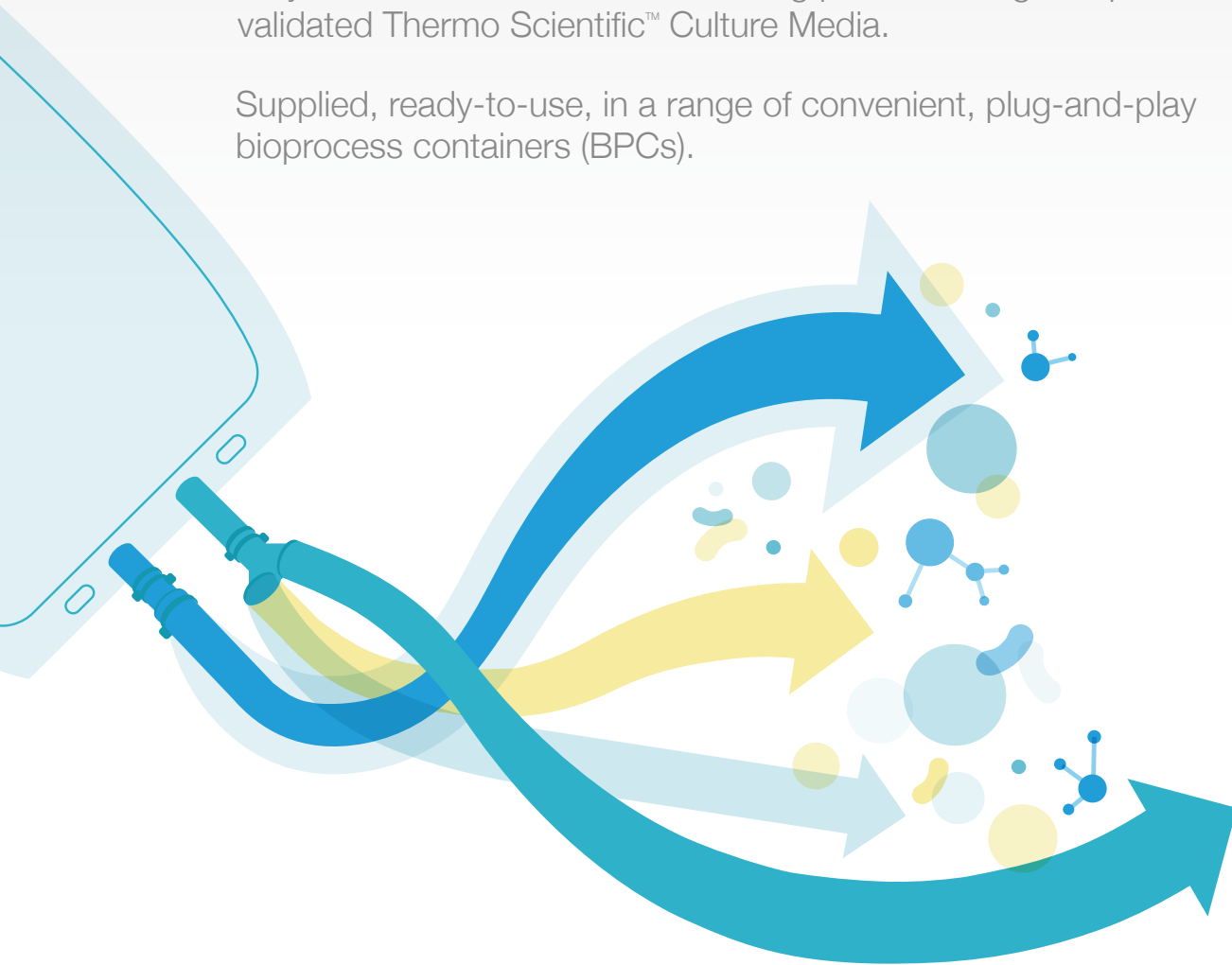

Complete confidence in 1, 2, 3

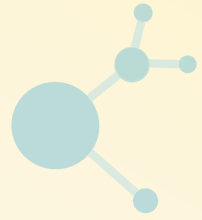

- 1** Achieve effective process simulations, with minimized risk of contamination
 - Innovative BPC design and connector options for easy process integration
- 2** Minimize preparation time
 - Standard, ready-to-use format
 - No water for injection
 - No component weighing
 - No mixing
 - No filter blockages
 - No vessel cleaning
- 3** Satisfy regulatory compliance
 - Pharmacopoeial-compliant Tryptone Soya Broth
 - Animal Derived Component Free (ADCF) vegetable peptone-based formulations available
 - No Mycoplasma contamination
 - Qualified microbial recovery performance
 - Comprehensive validation package for media, BPC, components and filling process



Media fills made easy


1 Achieve effective process simulation, with minimized risk of contamination

Using ready-prepared Thermo Scientific™ Cold Filterable Tryptone Soya Broth or Cold Filterable Vegetable Peptone Broth, supplied in Thermo Scientific Bioprocess Containers, allows you to :

- Evaluate the aseptic assembly and operation of critical, sterile equipment
- Qualify the operator
- Demonstrate that the environmental controls meet basic, sterile aseptic processing requirements

Efficient process simulations ensure the effectiveness of the filling process and minimize production downtime.

The innovative bag and connector design allow you to plug-and-play, efficiently replicating your filling processes while eliminating the risk of contamination.

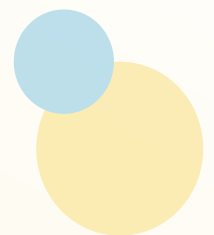



2 Minimize preparation time

Contaminated media may lead to false-positive results, leading to delays in reinstating production and additional costs.

Utilizing our experience as leaders in the manufacture of both dehydrated culture media and bioprocess containers, together with our liquid filling capabilities, we combine our cold filterable dehydrated media with water for injection (WFI) to provide you with an optimal placebo.

The ready-to-use media saves time and reduces labor costs by eliminating:

- Filter blockages
- Mixing
- Sterilization
- Vessel cleaning

The result is high-quality liquid media with a 12-month shelf life from date of manufacture.


**It's pure &
simple**

Convenient connections

- 1 ReadyMate™ aseptic connector may be used on its own or to make a sterile union to a tubing set to allow additional aseptic or steam connections



- 2 1½ inch sanitary flange connector

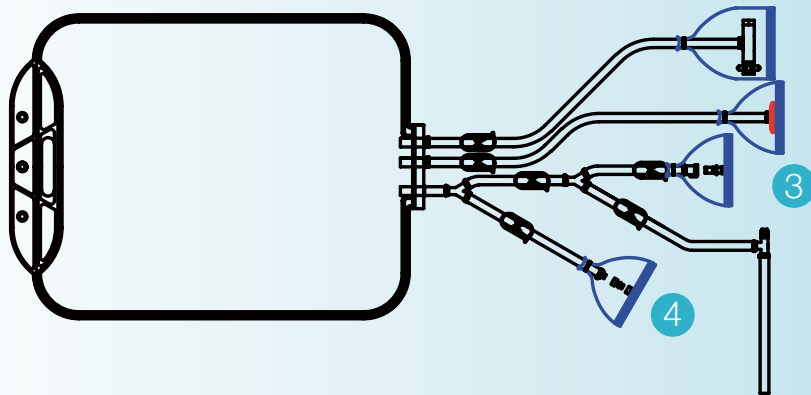


CHOOSE
READYMATE™
FOR EASY
CONNECTIONS
TO

- 3 CPC™ quick connect/disconnect system



- 4 SmartSite™ connector for needle-free additions, or sampling using a luer lock-style syringe



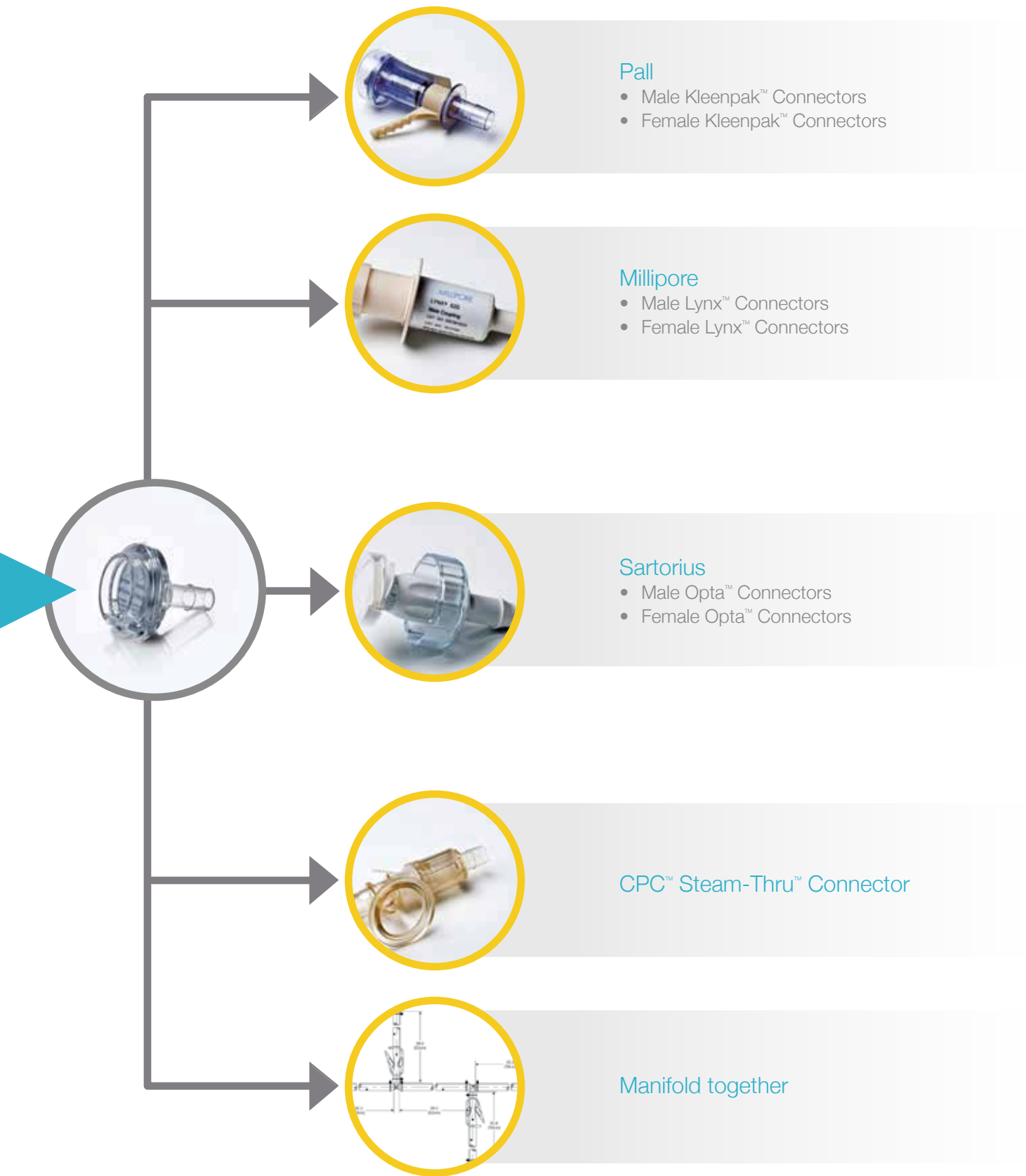

It's reassurance

you can trust

3 Satisfy regulatory requirements


Our gamma-irradiated, Cold Filterable Tryptone Soya Broth (TSB) and ADCF Vegetable Peptone Broth were formulated in our research and development laboratories, specifically for media fill trials.


Selected to maximize recovery and growth of microbial contaminants, ingredients used in our media are sourced for regulatory compliance, traceability and performance.

The media are gamma irradiated using a qualified cycle to ensure conformance to bioburden and Mycoplasma specifications.

Performance of our cold filterable TSB is tested according to specifications for the growth of control microorganisms laid down by the European, British, US and Japanese Pharmacopoeias.

Where an ADCF product is required, our vegetable peptone broth is a performance conformant alternative to Tryptone Soya Broth.

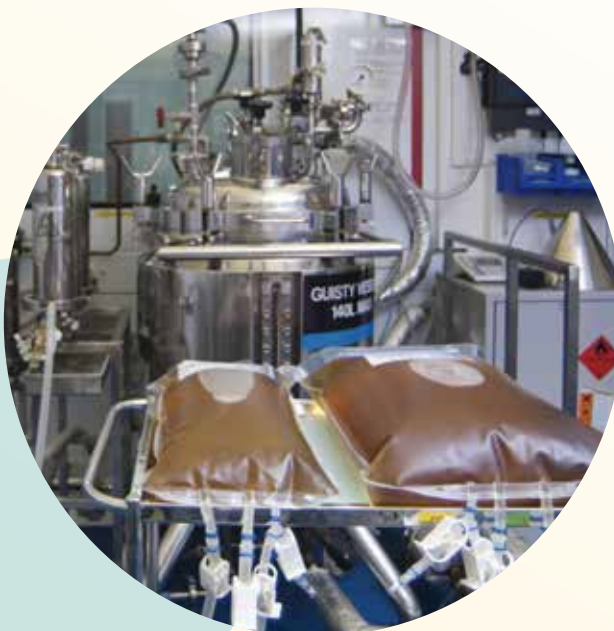

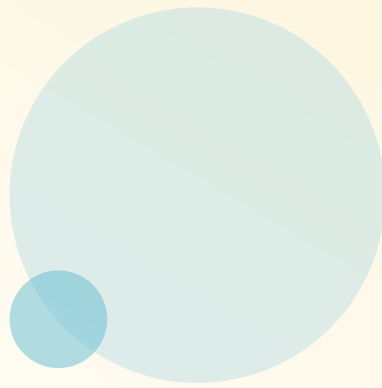

4 Save time & labor costs

Media fill trials must simulate accurately the filling production processes. Standard tryptone soya broths may cause filter blockages that can lead to interruptions in the process simulation.

Preparation, autoclaving and transporting of bulk media around the facility and time-consuming cleaning are now things of the past.

Our gamma-irradiated, cold-filterable media in 1-, 10- and 20-liter bioprocess containers are convenient and easy to use. They provide optimum filterability, with no need to heat the medium or change filters during the simulation, so that you save even more valuable time.

A comprehensive validation package covers our media, the BPC and its components, and the filling process.



Ordering information

1-liter options

Medium	Format	Order Code
Cold Filterable Tryptone Soya Broth	1L BPC	BP1065A
Cold Filterable Vegetable Peptone Broth	1L BPC	BP0104A

The above products are supplied with these connectors as standard :

- i) LUER lock female ¼"
- ii) SmartSite™ , needle-free valve additional sampling port

10- and 20-liter options

Medium	Format	Order Code
Cold Filterable Tryptone Soya Broth	10L BPC	BP1065C
	20L BPC	BP1065E
Cold Filterable Vegetable Peptone Broth	10L BPC	BP0104C
	20L BPC	BP0104E

The above products are supplied with these connectors as standard :

- i) SmartSite™ , needle-free valve additional sampling port
- ii) Quick connect, ¼" body (MPC)
- iii) 1½" triclamp
- iv) ReadyMate™ aseptic connector

Tubing Sets & Manifold Options

ReadyMate™ aseptic connectors may be used on their own, or to make a sterile union to a tubing set, to allow additional aseptic or steam connections.

The following tubing sets are available separately:

Product	Order Code
ReadyMate™ to male Lynx™	BP0010A
ReadyMate™ to female Lynx™	BP0020A
ReadyMate™ to male KPC	BP0030A
ReadyMate™ to female KPC	BP0040A
ReadyMate™ to male Opta	BP0050A
ReadyMate™ to female Opta	BP0060A
ReadyMate™ to Steam-Thru	BP0070A
ReadyMate™ to 2x ReadyMate™ “Y piece”	BP0080A
ReadyMate™ to 5x ReadyMate™ manifold	BP0090A



Thermo
SCIENTIFIC
A Thermo Fisher Scientific Brand

thermoscientific.com/BPCs

2015 © Thermo Fisher Scientific. All rights reserved. Kleenpak is a trademark or registered trademark of the Pall Corporation. Lynx is a trademark or registered trademark of Merck KGaA. Opta is a trademark or registered trademark of Sartorius Stedim Biotech. SmartSite is a trademark or registered trademark of Cardinal Healthcare Inc. ReadyMate is a trademark or registered trademark of GE Healthcare Companies. OPC and Steam-Thru are trademarks of Colder Products Company. All other trademarks are the property of Thermo Fisher Scientific and its subsidiaries.

Contact Information :
microbiology@thermofisher.com
USA +1 800 255 6730
International +44 (0) 1256 841144

993-090 Rev2
LT2018A Rev2
March 2015

代理商

啟新生物科技有限公司
www.cmp-micro.com

台北總公司 新北市新莊區新北產業園區五工五路21號 電話 (02) 2298-1823 傳真 (02) 2298-8100
台中分公司 台中市龍井區大同街241號 電話 (04) 2631-7419 傳真 (04) 2652-8275
高雄分公司 高雄市三民區建功路272巷6號 電話 (07) 392-1220 傳真 (07) 392-0748

C150529

Creative
CREATIVE LIFESCIENCES
啟新生技
唯有啟新 · 專業放心